# IRI

Immuno Enzymometric Assay

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**ST**  
E テスト「TOSOH」Ⅱ (IRI)

**ST AIA-PACK IRI**

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**en** Instructions For Use in multiple languages are available on our website: [www.tosohbioscience.eu](http://www.tosohbioscience.eu) under the "Service & Support" section. A paper version can be obtained by fax +32 13 66 47 49 or e-mail: [info.raqa.eu@tosoh.com](mailto:info.raqa.eu@tosoh.com).

⚠ This product contains human or animal source materials and should be treated as potentially infectious.

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### Attention

**For North and South American Customers: Please refer to the AIA-AAM Docs on CD for the appropriate information.**

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## ST AIA-PACK IRI

For Quantitative Measurement of insulin (IRI) in Serum, Heparinized Plasma or EDTA Plasma

### NAME AND INTENDED USE

ST AIA-PACK IRI is designed for IN VITRO DIAGNOSTIC USE ONLY for the quantitative measurement of insulin (IRI) in human serum, heparinized plasma or EDTA plasma on TOSOH AIA System Analyzers.

### SUMMARY AND EXPLANATION OF TEST

Insulin, the antidiabetic hormone, is produced in the pancreatic B cell as a large preproinsulin containing 109 amino acid residues with a molecular weight of approximately 11,500 (1,2). This peptide is rapidly converted by cleavage to proinsulin consisting of 86 amino acid residues with a molecular weight of approximately 9,000 and is stored within the B cell secretory granules (3). Equimolar quantities of insulin (51 amino acids with a molecular weight of approximately 6,000) and C-peptide (31 amino acids, molecular weight of approximately 3,000) are produced through proteolytic cleavage and then secreted along with a small amount of proinsulin (4). Insulin exists in polymeric forms depending on the pH and zinc content. The monomeric insulin molecule is composed of two polypeptide chains, "alpha" and "beta" which are connected by two interchain disulfide bridges of cystine (3).

Insulin is degraded in most tissues and has a plasma half-life of 7 - 15 minutes in man. It is also rapidly and completely inactivated in the gastrointestinal tract. Daily production of insulin in a healthy adult is 40 - 50 units. After binding to its specific receptors on target cell membranes, insulin acts as an anabolic and antidiabetic hormone, influencing the rates of carbohydrate, lipid, protein and electrolyte metabolism (5-7).

Insulin release is stimulated by glucose. A failure to respond to this glucose stimulus may be one of the fundamental defects in human diabetes (8-10). The factors stimulating or inhibiting insulin release or the factors decreasing the tissue response to insulin have been well documented.

### PRINCIPLE OF THE ASSAY

The ST AIA-PACK IRI is a two-site immunoassay performed entirely in the AIA-PACK test cups. Insulin present in the test sample is bound with monoclonal antibody immobilized on a magnetic solid phase and enzyme-labeled monoclonal antibody in the AIA-PACK test cups. The magnetic beads are washed to remove unbound enzyme-labeled monoclonal antibody and are then incubated with a fluorogenic substrate, 4-methylumbelliferyl phosphate (4MUP). The amount of enzyme-labeled monoclonal antibody that binds to the beads is directly proportional to the insulin concentration in the test sample. A standard curve is constructed, and unknown sample concentrations are calculated using this curve.

### MATERIAL PROVIDED (ST AIA-PACK IRI, Cat. No. 0025260)

Plastic test cups containing lyophilized eight magnetic beads coated with anti-insulin mouse monoclonal antibody and 100 µL of anti-insulin mouse monoclonal antibody (to human insulin) conjugated to bovine alkaline phosphatase with sodium azide as a preservative.

### MATERIALS REQUIRED BUT NOT PROVIDED

The following materials are not provided but are required to perform insulin analysis using the ST AIA-PACK IRI (Cat. No. 0025260) on the TOSOH AIA Systems Analyzers. They are available separately from TOSOH.

Materials	Cat. No.
AIA Nex-1A or AIA-21	0018539
AIA Nex-1A or AIA-21 LA	0018540
AIA-1800 ST	0019836
AIA-1800 LA	0019837
AIA-2000 ST	0022100
AIA-2000 LA	0022101
AIA-600 II	0019014
AIA-600 II BCR	0019328
AIA-360	0019945
AIA-PACK SUBSTRATE SET II	0020968
AIA-PACK SUBSTRATE REAGENT II	
AIA-PACK SUBSTRATE RECONSTITUENT II	
AIA-PACK IRI CALIBRATOR SET	0020360
CALIBRATOR (1) 0 µU/mL	
CALIBRATOR (2) 20 µU/mL (approx.)	
CALIBRATOR (3) 40 µU/mL (approx.)	
CALIBRATOR (4) 80 µU/mL (approx.)	
CALIBRATOR (5) 160 µU/mL (approx.)	
CALIBRATOR (6) 340 µU/mL (approx.)	
AIA-PACK IRI SAMPLE DILUTING SOLUTION	0020560
AIA-PACK WASH CONCENTRATE	0020955
AIA-PACK DILUENT CONCENTRATE	0020956
SAMPLE CUPS	0018581
AIA-PACK DETECTOR STANDARDIZATION TEST CUP	0020970
AIA-PACK SAMPLE TREATMENT CUP	0020971

Additional Requirements for AIA Nex-1A / AIA-21 only:

PIPETTE TIPS	0018552
PRELOADED PIPETTE TIPS	0018583
Additional Requirements for AIA-600 II, AIA-1800 and AIA-2000:	
PIPETTE TIPS	0019215
TIP RACK	0019216
PRELOADED PIPETTE TIPS	0022103

Only materials obtained from TOSOH should be used. Materials obtained elsewhere should not be substituted since assay performance is characterized based strictly on TOSOH materials.

### WARNINGS AND PRECAUTIONS

- The ST AIA-PACK IRI is intended for in vitro diagnostic use only.
- Test cups from different lots or different assays should not be mixed within a tray.
- The ST AIA-PACK IRI contains sodium azide, which may react with lead or copper plumbing to form potentially explosive metal azides. When disposing of such reagents, always flush with large volumes of water to prevent azide build-up.

- Human serum is not used in the preparation of this product; however, since human specimens will be used for samples and other quality control products in the lab may be derived from human serum, please use standard laboratory safety procedures in handling all specimens and controls.
- Do not use beyond the expiration date.
- The ST AIA-PACK IRI has been designed so that the high dose "hook effect" is not a problem for the vast majority of samples. Samples with insulin concentrations between 320 and 20,000 µU/mL will read > 320 µU/mL. The "hook effect" phenomenon may occur at insulin concentrations > 20,000 µU/mL.
- TOSOH AIA-1200 series or AIA-600 Immunoassay Analyzers can NOT be used to perform the ST AIA-PACK IRI assay.

### STORAGE AND STABILITY

All unopened materials are stable until the expiration date on the label when stored at the specified temperature.

Materials	Cat. No.
2-8°C:	
ST AIA-PACK IRI	0025260
AIA-PACK IRI CALIBRATOR SET	0020360
AIA-PACK IRI SAMPLE DILUTING SOLUTION	0020560
AIA-PACK SUBSTRATE SET II	0020968
AIA-PACK WASH CONCENTRATE	0020955
AIA-PACK DILUENT CONCENTRATE	0020956
1-30°C:	
AIA-PACK DETECTOR STANDARDIZATION TEST CUP	0020970
AIA-PACK SAMPLE TREATMENT CUP	0020971

ST AIA-PACK IRI test cups may be stored at 18-25°C for up to 1 day. Calibrators and Sample Diluting Solution should be used within 1 day and 7 days of opening or reconstituting, respectively, provided the vials are kept tightly sealed and refrigerated at 2-8°C. Reconstituted substrate solution is stable for 3 days at 18-25°C or 7 days at 2-8°C. Working diluent and wash solutions are stable for 30 days at 18-25°C. Reagents should not be used if they appear cloudy or discolored.

### SPECIMEN COLLECTION AND HANDLING

- Serum, heparinized plasma or EDTA plasma is required for the assay. Citrated plasma SHOULD NOT BE USED.
- If using serum, a venous blood sample is collected aseptically without additives. Store at 18-25°C until a clot has formed (usually 15-45 minutes), then centrifuge to obtain the serum specimen for assay.
- If using heparinized plasma or EDTA plasma, a venous blood sample is collected aseptically with designated additive. Centrifuge and separate plasma from the packed cells as soon as possible.
- Specimen types should not be used interchangeably during serial monitoring of an individual patient. Measured concentrations may vary slightly between sample types in certain patients.
- Samples may be stored at 2-8°C for up to 24 hours prior to analysis. If the analysis cannot be done within 24 hours, the sample should be stored frozen at -20°C or below for up to 60 days.
- Repeated freeze-thaw cycles should be avoided. Turbid serum samples or samples containing particulate matter should be centrifuged prior to testing. Prior to assay, bring frozen samples to 18-25°C slowly and mix gently.
- The sample required for analysis is 50 µL.

### PROCEDURE

For the AIA Nex-1A / AIA-21, AIA-600 II, AIA-1800, AIA-2000 and AIA-360, please refer to their Operator's Manual for detailed instructions.

#### I. Reagent Preparation

##### A) Substrate Solution

Bring all reagents to 18-25°C before preparing the working reagent. Add the entire contents of the AIA-PACK SUBSTRATE RECONSTITUENT II (100 mL) to the lyophilized AIA-PACK SUBSTRATE REAGENT II and mix thoroughly to dissolve the solid material.

##### B) Wash Solution

Add the entire contents of the AIA-PACK WASH CONCENTRATE (100 mL) to approximately 2.0 L of CAP Class I or NCCLS Type I Reagent Grade water, mix well, and adjust the final volume to 2.5 L.

##### C) Diluent

Add the entire contents of the AIA-PACK DILUENT CONCENTRATE (100 mL) to approximately 4.0 L of CAP Class I or NCCLS Type I Reagent Grade water, mix well, and adjust the final volume to 5.0 L.

#### II. Calibration Procedure


##### A) Calibration Curve

The calibrators for use with the ST AIA-PACK IRI have been standardized on WHO 1st IRP 66/304 (1974).

The calibration curve for ST AIA-PACK IRI is stable for up to 90 days. Calibration stability is monitored by quality control performance and is dependent on proper reagent handling and TOSOH AIA System maintenance according to the manufacturer's instructions.

Recalibration may be necessary more frequently if controls are out of the established range for this assay or when certain service procedures are performed (e.g. temperature adjustment, sampling mechanism changes, maintenance of the wash probe or detector lamp adjustment or change). For further information regarding instrument operation, consult the TOSOH AIA System Operator's Manual.

A sample calibration curve from AIA-1800 follows and shows the algorithm used for calculating results.
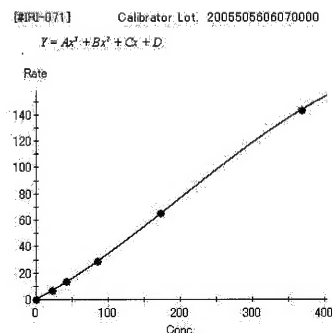
#### B) Calibration Procedure

1. Refer to the appropriate TOSOH AIA System Operator's Manual for the procedural instructions.
2. Verify that both the calibrator lot and concentration numbers have been correctly entered into the software.
3. CALIBRATOR (1) for ST AIA-PACK IRI is provided ready for use. TOSOH recommends that all calibrators should be run in triplicate.
4. CALIBRATOR (2)-(6) for ST AIA-PACK IRI are lyophilized. All levels should be reconstituted with 1.0 mL of CAP Class I or NCCLS Type I Reagent Grade water. TOSOH recommends that all calibrators should be run in triplicate.

#### C) Calibration Acceptability Criteria

1. The mean rate for the CALIBRATOR (1) should be  $< 3.0 \text{ nmol/(L}\cdot\text{s)}$ .
  2. Since there is a direct relationship between concentration and rate, the rate should increase as the concentration increases.
  3. The replicate values should be within a 10% range.
- D) Calibration Review and Acceptance
1. Review the calibration curve carefully, using the criteria listed above.
  2. Edit the calibration if necessary, then accept the calibration.

For further information regarding calibration, consult the TOSOH AIA System Operator's Manual.

#### III. Quality Control Procedure

##### A) Commercially Available Controls

Commercially available controls should be run at least once per day. It is recommended that at least two levels of controls, normal and abnormal, be used. Laboratory policy for this particular assay designates the following:

Control Material: \_\_\_\_\_  
Frequency: \_\_\_\_\_

Lot number of control material, acceptable limits, and corrective action to be taken if controls do not meet laboratory criteria will be found in a separate quality control document maintained by the laboratory.

##### B) Quality Control Procedure

1. Assay quality control specimens as instructed in the specific Operator's Manual for your analyzer. In addition, refer to the TOSOH AIA System Operator's Manual for detailed instructions on defining and editing the files.
2. Quality control material to be run with this assay is defined by individual laboratory policy.

#### IV. Specimen Processing

##### A) Preparation

Following specific instructions in the Operator's Manual for the analyzer, place samples on the instrument appropriately. Barcoded primary tubes as well as sample cups can be run on the AIA Nex-1A / AIA-21, AIA-600 II, AIA-1800, AIA-2000 and AIA-360.

##### B) Assay Procedure

1. Ensure a sufficient quantity of ST AIA-PACK IRI test cups for the number of samples to be run.
2. Load patient samples as instructed in the Operator's Manual and proceed with analysis. Note: The AIA Nex-1A / AIA-21, AIA-600 II, AIA-1800 and AIA-2000 will require AIA-PACK SAMPLE TREATMENT CUPS if onboard dilutions are utilized.

#### PROCEDURAL NOTES

1. Lyophilized substrate must be completely dissolved.
2. Ligand assays performed by the TOSOH AIA System Analyzers require that the laboratory use water designated by the College of American Pathologists as Class I or by NCCLS as Type I. Water should be tested at least once per month and should be free of particulate matter including bacteria. The pH of the water should also be routinely tested. For further information, consult the NCCLS document "Preparation and Testing of Reagent Water in the Clinical Laboratory," NCCLS Document C3-A3, Volume 11 No. 13, originally approved as a guideline by NCCLS in October 1997.
3. If a specimen insulin concentration is found to be greater than the upper limit of the assay range,  $320 \mu\text{U/mL}$ , the specimen should be diluted with the AIA-PACK IRI SAMPLE DILUTING SOLUTION and reassayed according to the Assay Procedure. The recommended dilution for specimens containing greater than  $320 \mu\text{U/mL}$  is 1:10 or 1:100. It is desirable to dilute the specimen so that the diluted specimen reads between 10 and  $300 \mu\text{U/mL}$ . The dilution factor should be entered into the software. For further information on the dilution of specimens, refer to the TOSOH AIA System Operator's Manual.
4. The TOSOH AIA System Analyzers can store two different calibration curves for each analyte at one time. Therefore, up to two different lots of ST AIA-PACK IRI test cups can be used during the same run.
5. If the assay specifications for this test are not ready in the system software, the specifications must be entered under test code **058**.

#### CALCULATION OF RESULTS

The TOSOH AIA System Analyzers perform all sample and reagent handling operations automatically. The TOSOH AIA System Analyzers read the rate of fluorescence produced by the reaction and automatically convert the rate to insulin concentration in  $\mu\text{U/mL}$ .

For samples requiring dilution, the AIA Nex-1A / AIA-21, AIA-600 II, AIA-1800 and AIA-2000 will automatically perform dilutions and calculate results if the dilution factors are entered into the software. Dilution factors may be entered into the Test File, or pre-defined dilution factors may be selected in Specimen Processing.

#### EVALUATION OF RESULTS

##### Quality Control

In order to monitor and evaluate the precision of the analytical performance, it is recommended that commercially available control samples should be assayed according to the local regulations.

The minimum recommendations for the frequency of running internal control material are:

- After calibration, three levels of the internal control are run in order to accept the calibration curve.
- The three levels of controls are repeated when certain service procedures are performed (e.g. temperature adjustment, sampling mechanism changes, maintenance of the wash probe or detector lamp adjustment or change).
- After daily maintenance, at least two levels of the control should be run in order to verify the overall performance of the TOSOH AIA System Analyzers.

If one or more control sample value(s) is out of the acceptable range, it is necessary to investigate the validity of the calibration curve before reporting patient results. Standard laboratory procedures should be followed in accordance with the strict regulatory agency under which the laboratory operates.

#### LIMITATIONS OF THE PROCEDURE

- For diagnostic purposes, the results obtained from this assay should be used in conjunction with other data (e.g. symptoms, results of other tests, clinical impressions, therapy, etc.).
- Using ST AIA-PACK IRI, the highest measurable concentration of insulin in specimens without dilution is  $320 \mu\text{U/mL}$ , and the lowest measurable concentration in specimens is  $0.5 \mu\text{U/mL}$  (assay sensitivity).
- Although the approximate value of the highest calibrator is  $340 \mu\text{U/mL}$ , the exact concentration may be slightly different. The assay specification, ASSAY RANGE HIGH, should be defined as the upper limit of the assay range,  $320 \mu\text{U/mL}$ .
- Hemolyzed samples should not be used for the ST AIA-PACK IRI assay because hemolysis may falsely bring down the value due to insulin-degrading enzyme in red blood cells.
- Lipemia has an insignificant effect on the assay except in the case of gross lipemia where spatial interference may occur.
- Specimens from patients taking medicines and/or medical treatment may show erroneous results.
- Samples containing fibrin may exhibit either falsely elevated or falsely decreased results.

- Specimens from patients who have received preparations of mouse monoclonal antibodies for diagnosis or therapy may contain human anti-mouse antibodies (HAMA). Such specimens may show falsely elevated values when tested for insulin.
- Specimens in the presence of autoantibody to insulin may show discrepant results.
- For a more complete understanding of the limitations of this procedure, please refer to the SPECIMEN COLLECTION AND HANDLING, WARNINGS AND PRECAUTIONS, STORAGE AND STABILITY, and PROCEDURAL NOTES sections in this insert sheet.

#### EXPECTED VALUES

Each laboratory should determine a reference interval which corresponds to the characteristics of the population being tested. As with all diagnostic procedures, clinical results must be interpreted with regard to concomitant medications administered to the patient (11).

#### Reference Ranges

The interval given here was determined in serum samples from 200 apparently healthy Asian individuals.

Reference Interval =  $1.1 - 17.0 \mu\text{U/mL}$  ( $1.1 - 17.0 \text{ mU/L}$ )

#### Conversion Factors

Insulin concentrations in this application are in units of  $\mu\text{U/mL}$ . Conversion to SI units of  $\text{mU/L}$  may be made using the following equation:  
 $\text{mU insulin/L} = \mu\text{U insulin/mL} \times 1.0$

#### PERFORMANCE CHARACTERISTICS

##### ACCURACY

- a. Recovery: Three serum pools were spiked with three different levels of insulin and assayed before and after spiking.

Sample	Initial Value ( $\mu\text{U/mL}$ )	Insulin Added ( $\mu\text{U/mL}$ )	Expected Value ( $\mu\text{U/mL}$ )	Measured Value ( $\mu\text{U/mL}$ )	Percent Recovery (%)
Serum A1	8.97	259.1	268.1	264.1	98.5
	8.97	129.6	138.5	142.7	103.0
	8.97	64.8	73.7	72.6	98.5
Serum B1	7.67	259.1	266.8	260.1	97.5
	7.67	129.6	137.2	138.4	100.9
	7.67	64.8	72.4	71.6	98.9
Serum C1	10.4	259.1	269.5	264.4	98.1
	10.4	129.6	140.0	141.7	101.3
	10.4	64.8	75.2	76.1	101.2

- b. Dilution: Three serum samples containing high concentrations of insulin were serially diluted with AIA-PACK IRI SAMPLE DILUTING SOLUTION and assayed.

Sample	Dilution Factor	Expected Value ( $\mu\text{U/mL}$ )	Measured Value ( $\mu\text{U/mL}$ )	Percent Recovery (%)
Serum A2	none		226.8	
	7.5/10	170.1	173.4	101.9
	5.0/10	113.4	116.7	102.9
	2.5/10	56.7	56.7	100.0
Serum B2	1.0/10	22.7	21.0	92.4
	none		265.2	
	7.5/10	198.9	207.2	104.2
	5.0/10	132.6	143.2	108.0
Serum C2	2.5/10	66.3	70.4	106.2
	1.0/10	26.5	25.1	94.6
	none		290.9	
	7.5/10	218.2	226.2	103.7
	5.0/10	145.5	156.0	107.3
	2.5/10	72.7	75.9	104.3
	1.0/10	29.1	27.6	94.9

##### PRECISION

- a. Within run precision was determined using three controls in a total of 20 runs. Within each run, one set of duplicates per control was assayed. The mean of each duplicate was used to obtain the pooled standard deviation (SD), which was then used to calculate the coefficient of variation (CV).

Sample	Mean ( $\mu\text{U/mL}$ )	Pooled SD ( $\mu\text{U/mL}$ )	CV (%)
Serum A3	12.3	0.277	2.3
Serum B3	92.2	1.324	1.4
Serum C3	164.4	2.436	1.5

- b. Total precision was determined by the duplicate assay of three controls in 20 separate runs. The means of each run were used to calculate the pooled standard deviation (SD) and coefficient of variation (CV).

Sample	Mean ( $\mu\text{U/mL}$ )	Pooled SD ( $\mu\text{U/mL}$ )	CV (%)
Serum A3	12.3	0.559	4.6
Serum B3	92.2	2.425	2.6
Serum C3	164.4	4.669	2.8

##### CORRELATION

The correlation between serum (x) and heparinized plasma (y) on ST AIA-PACK IRI was carried out using 91 patient specimens.

Slope	0.975
y-Intercept	-0.704
Correlation Coefficient	0.999
Range of Samples	1.2 - 281.9
Number of Samples	91

The correlation between heparinized plasma (x) and EDTA plasma (y) on ST AIA-PACK IRI was carried out using 60 patient specimens.

Slope	0.991
y-Intercept	0.719
Correlation Coefficient	0.998
Range of Samples	1.4 - 240.7
Number of Samples	60

##### SPECIFICITY

The following substances were tested for cross-reactivity. The cross-reactivity (%) is the percentage of the compound which will be identified as insulin. If these compounds are present in the specimen at the same concentration as insulin, the final result will be increased by these percentages.

Compound	Cross-reactivity (%)
Bovine insulin	72
Porcine insulin	100
Ovine insulin	75
Human C-peptide	N.D.
Human proinsulin	2
(N.D. : not detectable)	

#### SENSITIVITY

The minimal detectable concentration (MDC) of insulin is estimated to be 0.5  $\mu$ U/mL. The MDC is defined as the concentration of insulin which corresponds to the rate of fluorescence that is two standard deviations from the mean rate of fluorescence of 20 replicate determinations of a CALIBRATOR (1).

#### INTERFERENCE

Interference is defined, for the purposes of this study, as recovery outside of 10% of the known mean concentration of the specimen after the following substances are added to human specimens.

- Free bilirubin (up to 17 mg/dL) and conjugated bilirubin (up to 18 mg/dL) do not interfere with the assay.
  - Lipemia, as indicated by triglyceride concentration (up to 1,600 mg/dL), does not interfere with the assay.
  - Ascorbic acid (up to 20 mg/dL) does not interfere with the assay.
  - Protein, as indicated by human albumin concentration (up to 5.0 g/dL), does not interfere with the assay.
  - Heparin (up to 100 U/mL) does not interfere with the assay.
  - EDTA (up to 10 mg/mL) does not interfere with the assay.
  - Hemoglobin could interfere with the assay particularly at low range.
- Hemolyzed samples should not be used for the ST AIA-PACK IRI assay because hemolysis may falsely bring down the value due to insulin-degrading enzyme in red blood cells.

#### REFERENCES

1. Wrenshall, G. A., et al., The Story of Insulin. Indiana University Press. Bloomington (1962).
2. Lavy, P. L., The Pancreatic Beta Cell Structure and Function. New Eng. J. Med. 276, 187 (1967).
3. Steiner, D. F., Proinsulin and the Biosynthesis of Insulin. New Eng. J. Med. 283, 522 (1970).
4. Levine, R., Mechanisms of Insulin Secretion. New Eng. J. Med., 283, 522 (1970).
5. Rieser, P., Insulin, Membranes and Metabolism. Williams and Wilkins, Baltimore (1967).
6. Ashmore, J., Weber, G., Hormonal Control of Carbohydrate Metabolism in Liver in "Carbohydrate Metabolism and its Disorders" (ed. Dickens, Randle and Wholam) Academic Press, London, New York (1968).
7. Samols, E., Immunological Aspect of Insulin on the Nature and Treatment of Diabetes Mellitus. (ed. Leibel, Wrenshall) Excerpta Medica Foundation, Amsterdam (1965).
8. Fajans, S. S., Conn, J. W., The Early Recognition of Diabetes Mellitus. Ann. N. Y. Acad. Sci. 82:208 (1959).
9. Cerasi, E., Luft, R., Plasma Insulin Response to Glucose Infusion in Healthy Subjects and in Diabetes Mellitus. Acta Endocr. (kbh) 55:278-304 (1967).
10. Cerasi, E., et al., Decreased Sensitivity of the Pancreatic beta-Cells to Glucose in Prediabetic and Diabetic Subjects. Diabetes 21: 224-234 (1972).
11. Young, D., Effects of Drugs on Clinical Laboratory Tests. 3rd Edition, Washington, DC, American Association for Clinical Chemistry Press (1990).

#### TOSOH CORPORATION

Shiba-koen First Bldg.  
3-8-2 Shiba, Minato-ku, Tokyo 105-8623  
Japan  
Phone :+81-3-5427-5181  
Fax :+81-3-5427-5220

#### TOSOH EUROPE N.V.

Transportstraat 4  
B-3980 Tessenderlo, Belgium  
Phone :+32-13-66 88 30  
Fax :+32-13-66 47 49

#### TOSOH BIOSCIENCE, INC.

6000 Shoreline Ct., Suite 101  
South San Francisco, CA94080, USA  
Phone :+1-650-615-4970  
Fax :+1-650-615-0415  
Phone : (800)248-6764  
Fax : (800)685-7595